Medaysis

Enable Innovation

Mouse Anti-CD10 [MD182]: MC0272, MC0272RTU7

Intended Use: For Research Use Only

Description: The common acute lymphoblastic leukemia antigen (CALLA/CD10) is a single-pass type II transmembrane metallo-endopeptidase that cleaves and inactivates a variety of peptide growth factors important for signal transduction including the enkephalins, bombesin and substance P. CD10 is expressed by a number of hematopoietic cells such as immature T and B cells, B cells of the germinal centers of lymphoid follicles and granulocytes. It also reacts with a variety of non-hematopoietic cells, including epithelial cells in GI tract and kidney tubular. In liver, the bile canaliculi shows a moderate to strong staining. CD10 has been used for the identification and classification of certain types of malignant lymphoma and leukemia. CD10 is expressed in a high percentage of cases of acute lymphoblastic leukemia, follicular lymphoma, Burkitt's lymphoma, some hematopoietic tumors, and chronic myelogenous leukemias in lymphoid blast crisis. It is also known to be a marker of endometrial stromal cells. It is helpful in differentiating endometrial stromal sarcoma (ESS) from uterine cellular leiomyoma (UCL) and uterine leiomyosarcoma (ULS).

Specifications:	
Clone:	MD182
Source:	Mouse
Isotype:	IgG2b/k
Reactivity:	Human
Immunogen:	Recombinant human CD10 protein fragment aa 54-181
Localization:	Membrane, cytoplasm
Formulation:	Purified antibody in PBS pH7.4, containing BSA and $\leq 0.09\%$ sodium azide (NaN3)
Storage:	Store at 2°- 8°C
Applications:	IHC


Applications: Package:

скиде		Catalog No.	Size
	CD10 Concentrated	MC0272	1 ml
	CD10 Prediluted	MC0272RTU7	7 ml

IHC Procedure*:

Positive Control Tissue:	Kidney, prostate
Concentrated Dilution:	50-200
Pretreatment:	Tris EDTA pH9.0, 15 minutes using Pressure Cooker, or 30-60 minutes using water bath at 95°-99°C
Incubation Time and Temp:	30-60 minutes @ RT
Detection:	Refer to the detection system manual

* Result should be confirmed by an established diagnostic procedure.

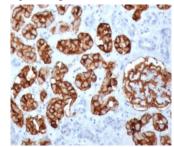

FFPE human kidney stained with anti-CD10 using DAB

References:

- 1. Surgical monotherapy may be a suitable therapeutic strategy for advanced collecting (Bellini) duct carcinoma: A case report and literature review. Zhu L, et al. Exp Ther Med 12:1181-1184, 2016.
- 2. Analysis of the clonal architecture of the human small intestinal epithelium establishes a common stem cell for all lineages and reveals a mechanism for the fixation and spread of mutations. Gutierrez-Gonzalez L, et al. J Pathol 217:489-96, 2009.
- 3. Changes in expression pattern of selected endometrial proteins following mesenchymal stem cells infusion in mares with endometrosis. Mambelli LI, et al. PLoS One 9:e97889, 2014.

Doc. 100-MC0272 Rev. C

Orders: <u>customercare@medaysis.com</u> Support: <u>techsupport@medaysis.com</u> Tel: 510-509-3153 <u>www.medaysis.com</u> © Medaysis Company